

## 6 BACTER FIX AND COOL

(Duration: ~35 min)

6.1 Notify MCC (Start of step)

Record MET: \_\_\_\_\_ / \_\_\_\_\_ : \_\_\_\_\_

CMP 6.2 ✓ Air intakes free of debris

GB 6.3 Perform GB Installation per Cue Card

CMP 6.4 CENTRIF 2 – OFF  
✓ CENTRIF 3 – OFF

### **NOTE**

Perform shaking and fixation in numerical order

Incubator 6.5 Transfer 4 Containers to GB (Shake vigorously 20x during transfer)

Cent 2      **BACTER 1**      to GB  
201              1–g

203              **BACTER 2**      to GB  
                    1–g

202              **BACTER 3**      to GB  
                    1–g

204              **BACTER 4**      to GB  
                    1–g

Incubator 6.6 Close, latch door

MF43H      6.7 Unstow:  
GB Tray      BACTER 9 (Large Screwdriver)  
                    BACTER 10 (6–Prong Plunger), transfer to GB

- GB 6.8 Perform GB Preparation per Cue Card  
(GB contents: BACTER 1, 2, 3, 4, BACTER 9 (Screwdriver), BACTER 10  
(6-Prong Plunger))  
Notify MCC start of GB Ops

GB Fix 8 Cassettes

6.9 **BACTER 1**  
1-g

Remove Cover  
Turn **red screw ccw** to open position "O"  
using Screwdriver

Push all **red** pistons (six) using 6-Prong  
Plunger

✓ All pistons down and red sleeve flush

Holding 6-Prong Plunger in place,  
turn **red screw cw** back to closed  
position

6.10 **BACTER 1**  
1-g

Return Cover and close


6.11 **BACTER 2** Repeat 6.9 – 6.10  
1-g

**NOTE**

There is no hard stop at the "O" position,  
turn **red screw ccw** no further than indicated on the red  
label (step 6.12)

6.12 **BACTER 3** Repeat 6.9 – 6.10  
1-g

6.13 **BACTER 4** Repeat 6.9 – 6.10  
1-g

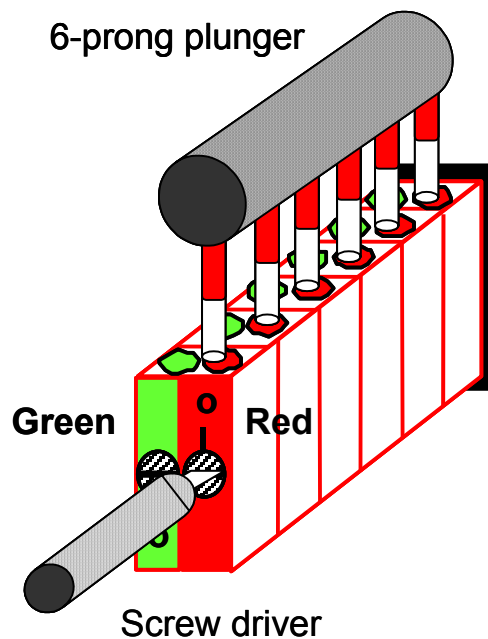

TOX LEVEL 2

- GB 6.14 GB sw – OFF  
Remove Gloving, temp stow  
Notify MCC end of GB Ops
- 6.15 Stow BACTER 9 (Screwdriver) and BACTER 10 (6-Prong Plunger) in GB Tray

|

**NOTE**

Perform shaking in numerical order

GB 6.16 Shake vigorously 20x **BACTER 1 to 4** and temp stow

Temp stow 6.17 Transfer 4 Containers to PTCU 2 (MF71H)

**CAUTION**

PTCU open minimal time for thermal integrity

**BACTER 1** to PTCU 2  
1-g

**BACTER 2** to PTCU 2  
1-g

**BACTER 3** to PTCU 2  
1-g

**BACTER 4** to PTCU 2  
1-g

GB 6.18 Perform GB Stow per Cue Card

6.19 Notify MCC (End of step)

Record MET: \_\_\_\_\_/\_\_\_\_\_:\_\_\_\_\_

7 INCUBATION 24 Hours (No crew required)

## 8 BACTER FIX AND COOL

(Duration: ~20 min)

8.1 Notify MCC (Start of step)

Record MET: \_\_\_\_\_ / \_\_\_\_\_ : \_\_\_\_\_

CMP 8.2 ✓ Air intakes free of debris

GB 8.3 Perform GB Installation per Cue Card

CMP 8.4 ✓CENTRIF 2 – OFF  
✓ CENTRIF 3 – OFF

### **NOTE**

Perform shaking and fixation in numerical order

Incubator 8.5 Transfer 4 Containers to GB (Shake vigorously 20x during transfer)

Rack        **BACTER 5**        to GB  
405        0-g

501        **BACTER 6**        to GB  
            0-g

Cent 3       **BACTER 7**        to GB  
301        0-g

302        **BACTER 8**        to GB  
            0-g

CMP 8.6 ✓CENTRIF 2 – OFF  
            ✓ CENTRIF 3 – OFF

Incubator 8.7 Close, latch door

MF43H       8.8 Unstow:  
GB Tray       BACTER 9 (Large Screwdriver)  
                 BACTER 10 (6-Prong Plunger), transfer to GB

- GB 8.9 Perform GB Preparation per Cue Card  
(GB contents: BACTER 5, 6, 7, 8, BACTER 9 (Screwdriver), BACTER 10 (6-Prong Plunger))  
Notify MCC start of GB Ops

GB Fix 8 Cassettes

8.10 **BACTER 5**  
0-g

Remove Cover  
Turn **red screw ccw** to open position "O"  
using Screwdriver

Push all **red** pistons (six) using 6-Prong  
Plunger

✓ All pistons down and red sleeve flush

Holding 6-Prong Plunger in place,  
turn **red screw cw** back to closed  
position


8.11 **BACTER 5**  
0-g

Return Cover and close

8.12 **BACTER 6** Repeat 8.10 – 8.11  
0-g

8.13 **BACTER 7** Repeat 8.10 – 8.11  
0-g

8.14 **BACTER 8** Repeat 8.10 – 8.11  
0-g

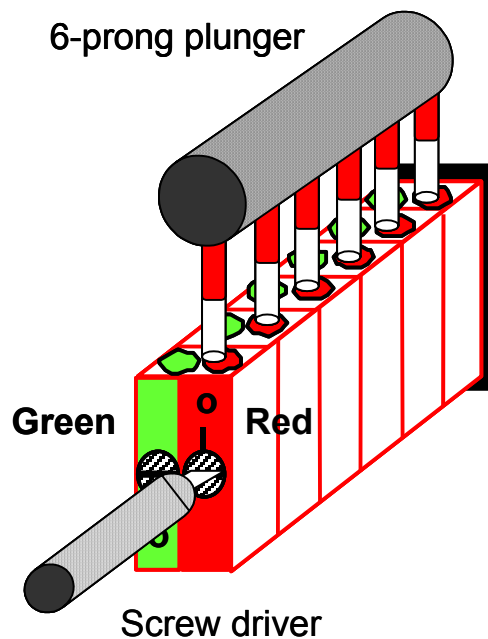

TOX LEVEL 2

- GB 8.15 GB sw – OFF  
Remove Gloving, temp stow  
Notify MCC end of GB Ops
- 8.16 Stow BACTER 9 (Screwdriver) and BACTER 10 (6-Prong Plunger) in GB Tray

|

**NOTE**

Perform shaking in numerical order

GB 8.17 Shake vigorously 20x **BACTER 5 to 8** and temp stow

Temp stow 8.18 Transfer 4 Containers to PTCU 2 (MF71H)

**CAUTION**

PTCU open minimal time for thermal integrity

**BACTER 5** to PTCU 2  
0-g

**BACTER 6** to PTCU 2  
0-g

**BACTER 7** to PTCU 2  
0-g

**BACTER 8** to PTCU 2  
0-g

GB 8.19 Perform GB Stow per Cue Card

8.20 Notify MCC (End of step)

Record MET: \_\_\_\_\_/\_\_\_\_\_:\_\_\_\_\_